
TERMOSIGILLATRICE

THERMOSEALER

THERMOSOUDEUSE

THERMOVERSIEGELUNGSGERÄT

TERMOSELLADORA

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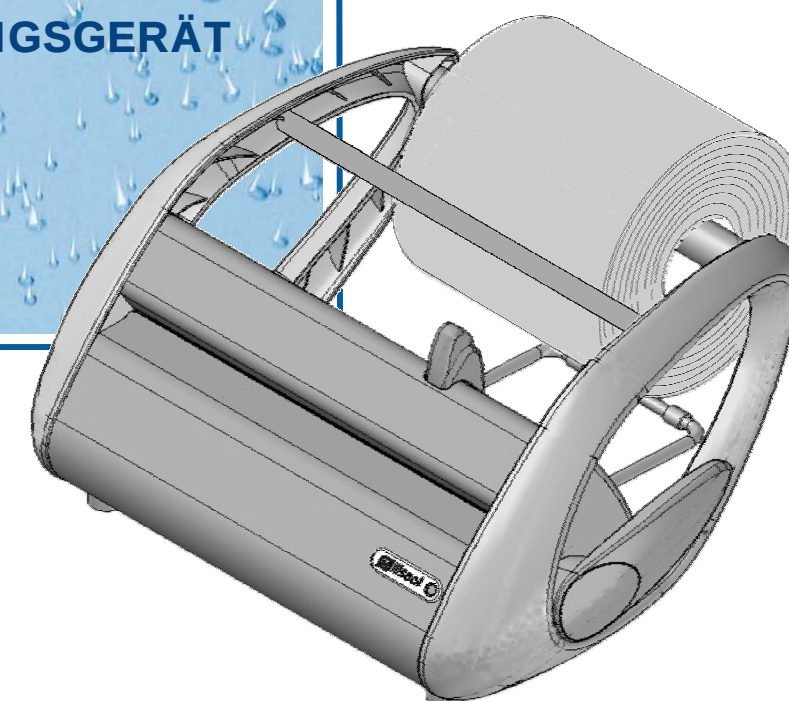THERMOSEALER OPERATING MANUAL

REVISIONS

The following table lists subsequent editions/revisions of the manual.
The "Description" field brief explains the subject of the latest revision.

Edition	Rev.	Date	Description
1	0	09-2005	First issue
1	1	06-2007	Added a paper-support rod

TABLE OF CONTENTS

INTRODUCTION	1
APPLICABLE EUROPEAN DIRECTIVES	1
INTENDED USE	1
PURPOSE OF THE MANUAL	2
GENERAL WARNINGS	2
CONTENTS OF THE PACKAGE	3
DIMENSIONS AND WEIGHT	3
DESCRIPTION OF THE ITEMS	3
PRODUCT PRESENTATION	4
MAIN PERFORMANCES	4
VIEWS	4
INSTALLATION	5
TECHNICAL DATA	5
GENERAL PRECAUTIONS	5
ELECTRICAL CONNECTION	6
MOUNTING THE ROLL	6
MOUNTING THE PAPER-SUPPORT ROD	6
OPERATION	7
PREPARING THE MATERIAL TO BE SEALED	8
PREMISE	8
MOVING THE MATERIAL	8
CLEANING THE TOOLS	8
WRAPPED MATERIAL	9
APPENDIX A – TECHNICAL CHARACTERISTICS	10
SUMMARY TABLE	10
SAFETY DEVICES	10

APPENDIX C - MAINTENANCE	11
GENERALS.....	11
CLEANING	11
BLADES REPLACEMENT	11
FUSES REPLACEMENT	11
APPENDIX D – TROUBLESHOOTING	12
APPENDIX F – ELECTRICAL DIAGRAM	13
APPENDIX G – DECLARATION OF CONFORMITY	14
APPENDIX Z - CUSTOMER SUPPORT	15

INTRODUCTION

Note, warning and danger symbols used through the manual

APPLICABLE EUROPEAN DIRECTIVES

INTENDED USE

Important notes

Dear Customer

Thank you for choosing a product from M.O.COM. Srl. We hope that you will find it completely satisfactory.

This manual describes all procedures for the correct installation and use of the equipment.

Should you have any questions or suggestions concerning this product or our support, do not hesitate to contact us.

NOTA

PAY SPECIAL ATTENTION; NOTES FOR PROPER INSTALLATION AND MAINTENANCE PROCEDURES AND USE METHODS.

WARNING

This symbol indicates a potential danger of injury. Follow the procedures described in the manual to avoid injuring the user and/or others.

DANGER

THIS SYMBOL SIGNALS A POSSIBLE DANGER FOR PEOPLE. FOLLOW THE SUGGESTED PROCEDURE IN ORDER TO PREVENT POSSIBLE DAMAGE TO USERS AND/OR OTHERS.

ENGLISH

The product described in this manual is made according to the highest safety standards and does not present any danger for the operator provided it is used according to the instructions reported below.

The product complies with the directives **89/336/CEE** and **73/23/CEE** (and following updates).

The unit is meant for professional use and can only be employed for the purpose it was designed for, that is the thermo-sealing and cutting of **sterilization paper/polypropylene rolls**.

NOTE

THE UNIT CANNOT WORK WITH ROLLS OF A MATERIAL OTHER THAN THE ONES STATED ABOVE.

WARNING

THE DEVICE MUST ONLY BE USED BY QUALIFIED PERSONNEL. IT MAY NOT BE USED OR HANDLED BY INEXPERT AND/OR UNAUTHORIZED PERSONNEL FOR ANY REASON.

NOTES

THE INFORMATION IN THIS MANUAL IS SUBJECT TO CHANGE WITHOUT NOTICE.

M.O.COM. SRL SHALL NOT BE RESPONSIBLE FOR DIRECT, INDIRECT, ACCIDENTAL OR CONSEQUENTIAL DAMAGE OR OTHER DAMAGES RELATED TO THE SUPPLY/USE OF THIS INFORMATION.

THIS DOCUMENT MAY NOT BE REPRODUCED, ADAPTED OR TRANSLATED, IN WHOLE OR IN PART, WITHOUT THE PRIOR, WRITTEN AUTHORIZATION OF M.O.COM. SRL

 IS REGISTERED TRADEMARKS OF M.O.COM. SRL

PURPOSE OF THE MANUAL

The manual purpose is to provide instructions for:

- General knowledge of the product;
- correct installation;
- its safe and efficient use;
- the treatment of tools and/or materials to submit to sterilization.

More the appendix includes:

- The technical features of the product;
- Maintenance operations;
- Troubleshooting and solutions;
- Other documentation

GENERAL WARNINGS

The product must **always** be used according to the procedures included in this manual and never for purposes other than the ones indicated.

WARNING

THE USER IS RESPONSIBLE FOR ALL LEGAL REQUIREMENTS RELATED TO THE INSTALLATION AND USE OF THIS PRODUCT. THE MANUFACTURER WILL NOT BE RESPONSIBLE FOR ANY BREAKAGE, MALFUNCTIONS, PROPERTY DAMAGE OR INJURY IN THE EVENT THAT THE PRODUCT IS NOT INSTALLED OR USED CORRECTLY.


Please observe the following precautions in order to avoid injury or property damage; in particular:

- Use **ONLY** recommended cleaning detergents.
- **Do not** pour water or other liquids or inflammable products on the equipment;
- **Do not** use the equipment in the presence of gas or explosive or inflammable vapors;

DANGER

BEFORE PERFORMING ANY MAINTENANCE OR CLEANING, DISCONNECT THE PLUG OF THE POWER SUPPLY CABLE FROM THE SOCKET.

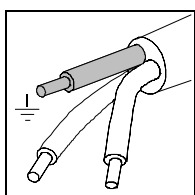- Make sure the electrical system is **grounded** conforming to current laws and/or standards;
- **Do not** remove any label or nameplate from the device; request new ones, if necessary.
- Use **only** original replacement parts.

WARNING



THE FAILURE TO OBSERVE THE ABOVE, RELEASES THE MANUFACTURER FROM ALL LIABILITY.

CONTENTS OF THE PACKAGE


DIMENSIONS AND WEIGHT

Package:

P Total weight about 8 Kg

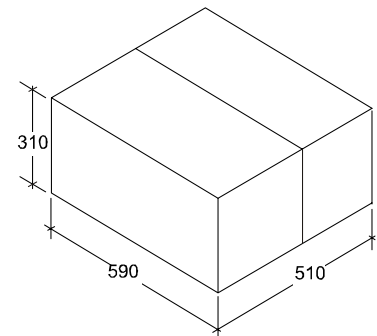
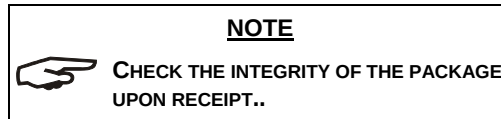
- Once the package is opened, check that:
- the supply matches the specifications of the order (see the accompanying document);
 - that there is no obvious product damage.

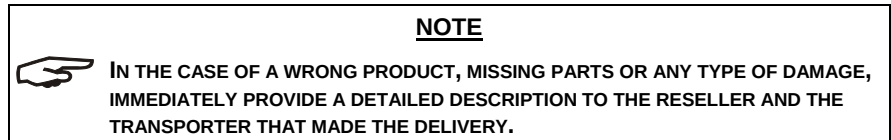ENGLISH

DESCRIPTION OF THE ITEMS

- * Thermosealer **Millseal**
- * Roll support (Ref. 1)
- * Paper-support rod (Ref. 2)
- * Power supply cable (Ref. 3)
- * Operating Manual (Ref. 4)
- * Warranty Certificate (Ref. 5)

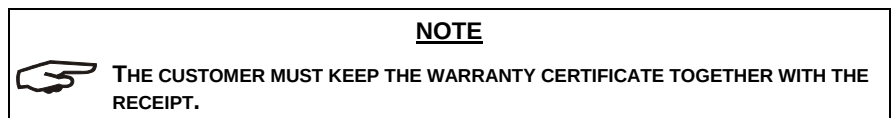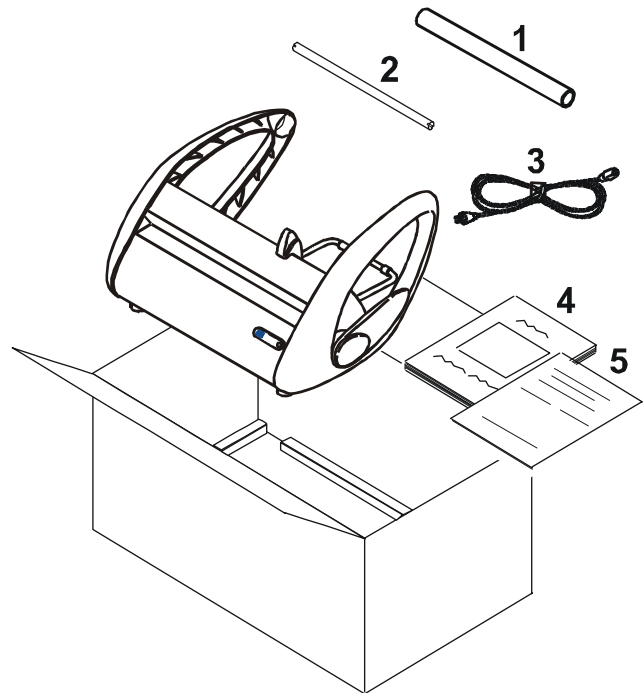PRODUCT PRESENTATION

Millseal is the new thermo-sealer by M.O.COM., which is positioned as a natural complement of the steam sterilizers range **Millennium** for professional use.

It is a unit with a pleasant design with a reduced footprint that combines ease of use and high reliability, also in case of intense usage.

Millseal represents therefore an interesting option for all users that need professional equipment with a good price/performance ratio.

MAIN PERFORMANCES

Millseal is an electronically controlled thermo-sealer with a sealing band of 12 mm fit to resist the sterilization processes of the modern sterilizers featuring fractioned pre-vacuum, and to assure in time the sterility of the packages.

Thanks to the temperature sensor integrated with the resistance, an effective heating control in the sealing area is obtained, ensuring high and stable performances.

The plastic and steel robust structure together with the ergonomic and elegant shape, allow to maintain a high standard of hygiene, reduced weight and size, as well as an exceptional strength.


All this has the precise goal of providing a product suitable for daily needs in a medical environment, also keeping into account the refined character of many surgeries interior design.

NOTE

FOR THE DESCRIPTION OF THE SECURITY DEVICES, REFER TO APPENDIX A (TECHNICAL CHARACTERISTICS)

VIEWS

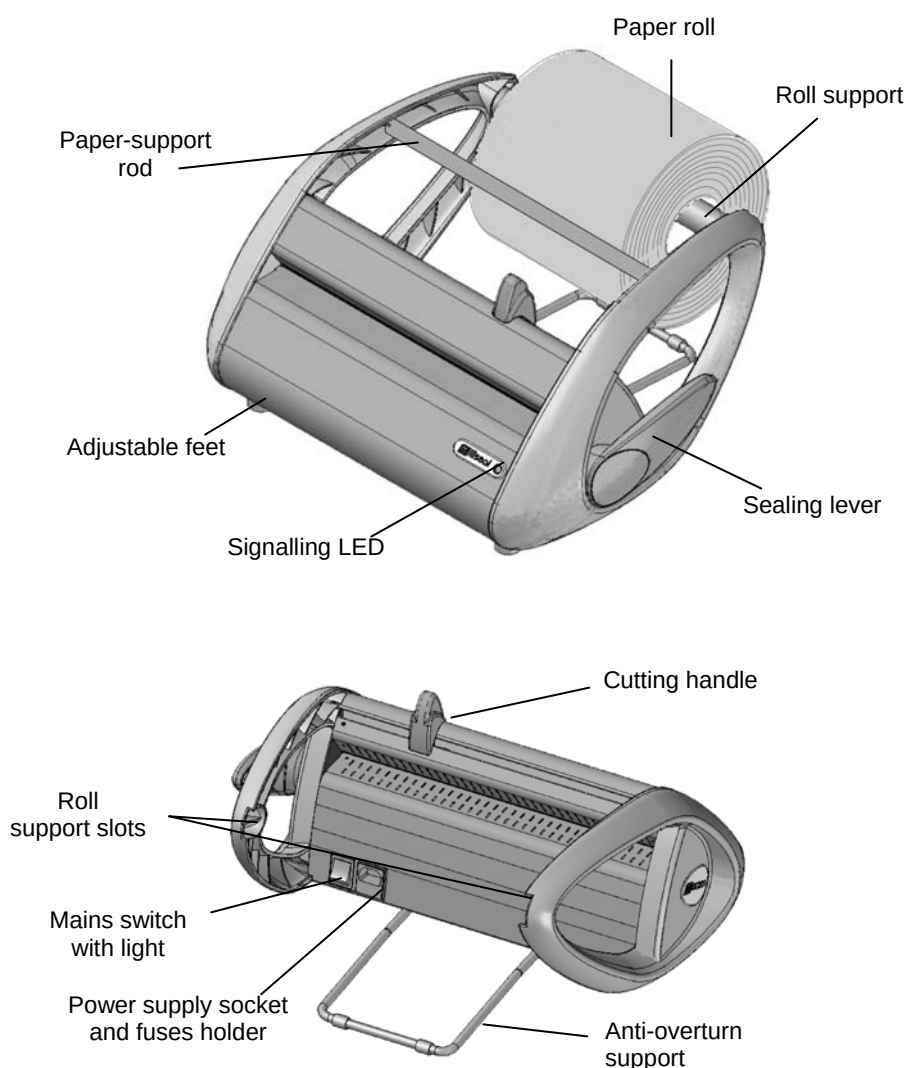INSTALLATION

For a good operation of the thermo-sealer the first and essential task is its correct installation; in this way possible malfunctions or damages to the unit will be avoided.

Therefore please read carefully the notes and warnings reported in the rest of this chapter.

NOTE

THE M.O.COM. TECHNICAL SERVICE (SEE APPENDIX Z) IS AT YOUR DISPOSAL, TO CLEAR ANY QUESTIONS OR PROVIDE FURTHER INFORMATION.

TECHNICAL DATA

Size and weight of the unit

- Max width 497 mm
- Max depth 394 mm (without roll)
- Max height 239 mm (without roll)
- Overall weight 7 kg (without roll)

Maximum dimensions of the paper/
polypropylene film roll Ø= 200, H = 300 mm

Power supply

Voltage: 220-240 Vac, 50 Hz

Power: 100 W – 0,45 A

The data plate is placed back on the unit.

ATTENTION

THE ELECTRICAL INSTALLATION MUST BE PROVIDED WITH EARTH CONNECTION.

GENERAL PRECAUTIONS

- Install the thermo-sealer on a flat surface, leaving enough room for its operation and adequate ventilation.
Makes sure that the resting plane is strong enough to support the unit and roll weight (max 12 kg);
- Do not install near sinks to avoid contact with water or liquids that could cause short circuits or situations potentially dangerous for the operator;
- Do not install in humid or poorly ventilated environments;
- Do not install in environment with the presence of inflammable and/or explosive **gas** or **steams**;
- Install the unit so that the power supply cable is not bent or pressed along its route up to the mains socket.

ELECTRICAL CONNECTION

In compliance with the standing regulations, the unit must be connected to a mains socket of the electrical installation with power rating suitable for the unit consumption and provided with earth connection.

The socket must be suitably protected upstream by a terminal board provided with electro thermal magnetic switch and differential switch.

Connect the supplied power cable to the mains socket on the rear of the equipment.

NOTE


ONLY USE THE POWER SUPPLY CABLE PROVIDED WITH THE UNIT; POSSIBLE PLUG REPLACEMENT IS UNDER THE FULL RESPONSIBILITY OF THE USER.

DIRECTLY CONNECT THE POWER SUPPLY CABLE TO THE SOCKET. Do NOT USE EXTENSIONS, ADAPTERS OR OTHER ACCESSORIES.


MOUNTING THE ROLL

Insert the paper-roll on the roll support and mount the assembly on the arm's roll support slots, and turned so that the paper unrolls downwards and with the polypropylene film upwards, as shown in the figure.

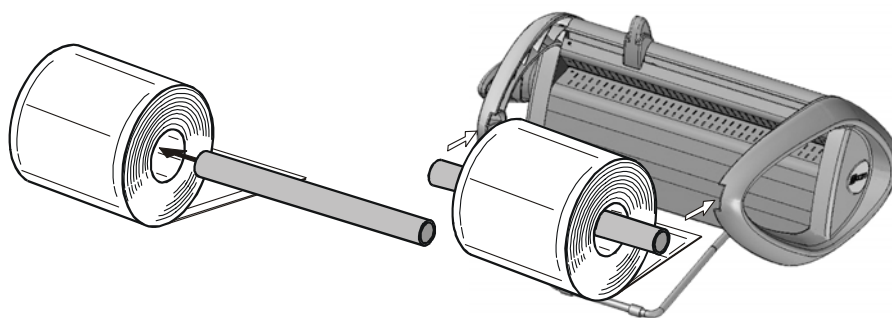MOUNTING THE PAPER-SUPPORT ROD

By practising a light bending, insert the paper-support rod between the two central and opposite ribs of the roll-support arms, as shown in the figure.

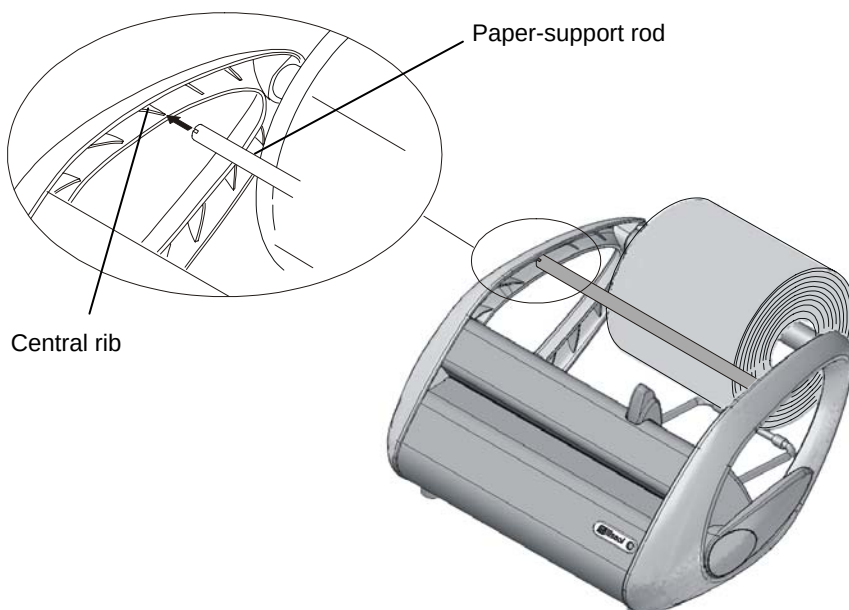OPERATION


Switch on the unit by means of the power supply switch placed in the unit rear. The switch illuminates while the front Led is still off.

The heating resistance is powered and the temperature gradually increases. When the temperature reaches the optimal value for a correct sealing, the front Led switches on (colour **GREEN**).

The thermo-sealing takes place in three steps:


1. INSERTING THE PAPER

Unroll the paper tape over the support rod, insert in the slot between the upper cutting carriage and the cover of the heating area, indicatively in the central zone, and extract it off the front for the desired length.

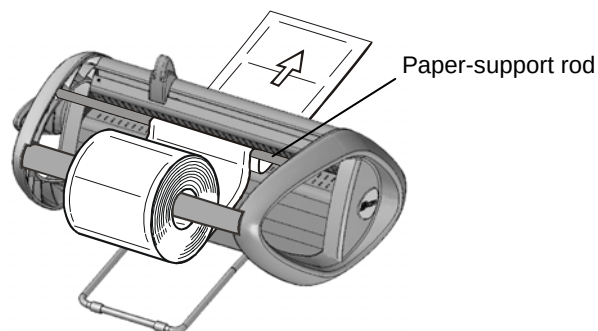2. SEALING ONE END AND CUTTING


Lift the sealing lever and keep it raised until you hear the beep and the Led switches colour (RED); always keeping the lever raised, move the cutting handle along the carriage. Return the lever to its original position. The LED goes back to **GREEN** and the acoustic signal stops. The envelope has been sealed at one end; take the envelope and proceed with the sealing of the other end.

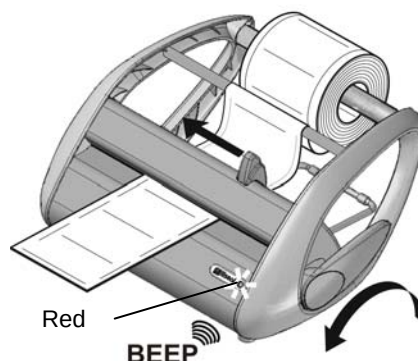3. PUTTING THE INSTRUMENT AND FINISHING THE ENVELOPE

Put the instruments in the envelope from the open side, following the instructions contained in the chapter "Preparing the material to be sealed" (see next page).

Insert the unsealed end of the envelope for a few centimetres from the front and, proceeding as described at point 2 above, carry out the second sealing and cut the unused paper part, if necessary.

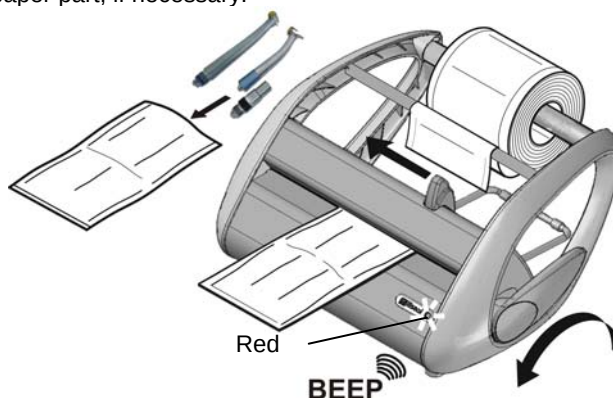PREPARING THE MATERIAL TO BE SEALED

PREMISE

The sterilization process can be considered as effective, reliable and repeatable on condition that the instruments and the materials in general are first properly treated and subsequently correctly positioned inside the sterilization chamber. This foresight is even more required when we consider **objects within an envelope**.

It must be stressed that organic residue or deposit of substances used in the medical practice inevitably become receptacle of micro-organisms and can hamper the steam sterilization process.

An improper arrangement of the load makes difficult and sometimes impossible the circulation of water steam on the material and/or inside its cavities, with the consequences that can be imagined.

Therefore we provide a few **basic** indications concerning these aspects, leaving the customers with the task of improving their knowledge on the matter in the best fashion.

MOVING THE MATERIAL

When **handling and moving** contaminated material, it is good practice to take the following precautions:

- Always wear rubber gloves with adequate thickness;
- Clean with a germicide detergent the hands protected with the gloves;
- Always use a tray for carrying the tools.
- Never carry the tools directly in your hands;
- Pay attention to any sharp or cutting parts to avoid the risk of infections;
- Separate the objects that must not undergo a sterilization or cannot stand the process;
- When you have finished handling the material, carefully wash your hands still gloved.

CLEANING THE TOOLS

All materials and/or tools to be sterilized must be perfectly cleaned and without trace of any kind of residue (organic and inorganic materials, paper fragments, cotton or gauze swabs, scale, etc.).

Otherwise, the sterilization process could turn out to be not effective, and the sterilizer and the tools themselves undergo irreversible damage.

For an effective cleaning:

1. After use, **immediately** rinse the tools with a jet of running water;
2. Divide the tools according to the type of material (carbon steel, stainless steel, brass, aluminium, chrome etc.) to avoid processes of electrolytic oxidation-reduction;
3. Carry out a washing with an ultrasonic cleaner containing a mix of water and germicide solution, carefully following the manufacturer's recommendations.
4. Carefully rinse and check that the cleaning is complete; if necessary repeat the washing cycle or clean manually.
5. Dry perfectly.

Before the sterilization, check the indications provided by the manufacturer of the tool/material, checking any incompatibility.

Meticulously follow the use procedures for detergent or disinfectant products and the instructions for the washing and automatic lubrication equipment.

WRAPPED MATERIAL

The use of sterilization paper envelopes allows both an optimal sterilization and to maintain in time the sterility condition.

Follow these general indications:

- Envelope the tools one by one or make sure that they are made of the same material;
- Close the envelope with the thermo-sealer.
- Do not use metal points, pins or similar since they do not warrant sterility;
- Arrange the envelopes so as to avoid the creation of air pockets that prevent the correct penetration and the removal of the steam.
- If possible, arrange the envelopes vertically using a suitable support.
- **Never** overlap the envelopes one over the other.

APPENDIX A – TECHNICAL CHARACTERISTICS

SUMMARY TABLE

Device	Thermo-sealer
Model	
Manufacturer	M.O.COM. S.r.l. Via delle Azalee, 1 20090 BUCCINASCO (MI) - ITALY
Power supply voltage	220/240 Vac
Frequency	50 Hz
Mains fuses (5 x 20 mm)	F 3,15A 250V
Nominal power	100 W
Insulation class	Class I
Installation category	Cat. II
Utilization	Internal use
Operating environment conditions	Temperature: +15°C ÷ +35°C Relative humidity: 80% (max), without condensation
External dimensions (without roll) (L x D x H)	497 x 394 x 239 mm
Net weight	7 kg
Type of sealing	With continuous cycle and constant temperature
Nominal operating temperature of the heating resistance	180°C
Depth of the sealing band	12 mm
Max sealing length	310 mm
Sealable material	paper /polypropylene film
Max diameter of the roll	200 mm
Orientation of the paper /polypropylene film	Polypropylene film oriented upward

SAFETY DEVICES

The thermo-sealer is equipped with the following safety devices, of which we list a short description along with their function:

- **Mains fuses** (see data in the above table)
Equipment protection from internal electrical faults.
Action: break the electrical power supply.
- **Temperature sensor for the heating resistance**
Protection against overheating
Action: interrupts the heating circuit.
- **Acoustic and visual indication for sealing completed**
Protection against envelope overheating.
Action: Led colour changeover (red-green) and acoustic signal for sealing completed.
- **Anti- finger crushing design**
Protection accident crushing of one's fingers during the lever operation, through fixed turret and mobile resistance built into the machine body.
- **Protected cutting knife**
Exclusive design with cutting knife totally inaccessible from outside.

APPENDIX C - MAINTENANCE

GENERALS

The equipment does not need any special care; it is only recommended that the user cleans it periodically.

In all cases carefully follow the procedure detailed below.

In case of replacement of equipment components, request and use ONLY ORIGINAL SPARE PARTS.

DANGER

BEFORE PERFORMING ANY MAINTENANCE OPERATION, DISCONNECT THE PLUG OF THE POWER SUPPLY CABLE FROM THE SOCKET.

ENGLISH

CLEANING

Clean the equipment externally with a clean cloth dampened with water and neuter detergent.

Perfectly dry the surfaces and remove any remaining resistant residue.

For the disinfection of external surfaces, always use weak detergents made up of a minimum percentage of sodium hypochlorite (or equivalent).


Always keep in mind the following **general warnings**:

- **Do not** submit the unit to direct water jets, under pressure or sprinkling. Any infiltration on the internal electrical components could irretrievably damage its operation;
- **Do not** use abrasive sponges, metal brushes (or other aggressive materials) or products for cleaning of metal, both solid and liquid;

BLADES REPLACEMENT

Periodically, or when you notice a difficulty in the sterilizer paper cutting, it is necessary to replace or rotate the position of the blades used for the envelope cutting.


- Remove the 4 fixing screws of the cutting carriage cover to access the paper cutting assembly.
- Remove the 3 nuts of plate "A" that keep into place the blades "B".
- Replace or invert the blade position.
- Reassemble the components following the inverse order of its disassembly.

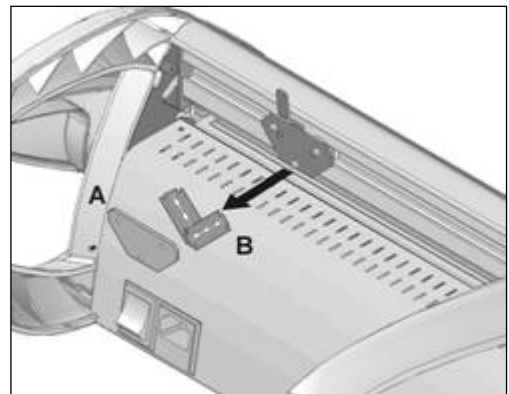FUSES REPLACEMENT

The fuses are housed in the fuses holder under the socket for the power supply cable.

- Remove the power supply cable from the socket .
- Extract the fuses drawer by a screwdriver.
- Replace the fuse with one of equal value.

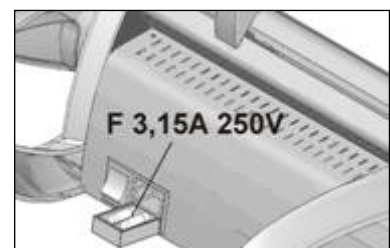APPENDIX D – TROUBLESHOOTING

If Your Thermo-sealer does **not** work correctly, carry out the following checks **before** contacting the Technical Service:

PROBLEM	POSSIBLE CAUSE	SUGGESTED SOLUTION
The thermo-sealer does not switch on.	The plug of the power supply cable is not inserted in the outlet.	Insert the plug.
	Lack of power at the power supply outlet.	Check the cause of the power failure and correct it.
	The magneto thermal switch and/or differential switch of the installation are OFF.	Put the switches to ON.
	Mains fuses broken	Replace the fuses (see <u>Appendix C - Maintenance</u>)
The Led does not switch on after a few minutes from the power on	The Led does not work.	Contact the Technical Support Service (see <u>Appendix Z</u>).
	Heating resistance broken.	
	Fault in the electronic card.	
The envelope sealing takes place but the Led remains green (instead of red).	Fault of the end-of-run microswitch.	Contact the Technical Support Service (see <u>Appendix Z</u>).
Unsatisfactory sealing of the envelope.	Insufficient pressure in the sealing area.	Contact the Technical Support Service (see <u>Appendix Z</u>).
	Silicon rubber profile damaged.	
Envelope cutting is difficult.	Consumed blades.	Replace or invert the position of blades (See <u>Appendix C - Maintenance</u>)

NOTE

IF THE PROBLEM PERSISTS, CONTACT THE TECHNICAL SUPPORT SERVICE (SEE **APPENDIX Z) REPORTING THE **THERMO-SEALER MODEL** AND ITS **SERIAL NUMBER**.**

THESE DATA CAN BE FOUND ON THE FACTORY PLATE PLACES ON THE UNIT REAR.

The diagram illustrates the electrical connections for the MCI 1024. At the bottom, the 'INPUT 220-240Vac' is shown with two main lines and a ground. Each line passes through a 3.15A fuse (labeled 'F 3.15A' and 'MAINS FUSES'). The ground line is connected to a common ground symbol. The two main lines then branch out: one goes to a terminal on the 'ELECTRONIC BOARD' (terminal 10), and the other goes to a terminal on the 'MCIRSWITCH' (terminal 1). The 'ELECTRONIC BOARD' is a large rectangle containing various components like capacitors, resistors, and a central square component. The 'MCIRSWITCH' is a vertical rectangle with a switch symbol. The 'HEATER' is represented by a zigzag line. Wires connect the board to the heater and the switch, and the switch to the heater. The board also has a terminal strip at the top with terminals 1 through 10.

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DECLARATION OF CONFORMITY

Application of the EEC Directives 93/42/CEE -89/336/CEE-73/23/CEE

Name of the Manufacturer : **M.O.COM. S.r.l.**
Manifattura Odontoiatrica Complementare

Address of the Manufacturer: **Via delle Azalee, 1 - 20090 Buccinasco (MI) – ITALIA**

Product description: **Thermo-sealer**

Model: 

Made in: **ITALY**

The undersigned hereby declares that the above mentioned goods entirely

conform

to the EEC Directives 89/336 - 73/23 (and following updates).

Reference Standards:	EN 55014-1	EN 55014-2
	EN 61000-3-2	EN 61000-3-3
	EN 61010-1	EN 61010-1-A2

15/09/2005

(Date)

Signature

Alfio VILLA
Name and Surname

Legal Agent
Function

APPENDIX Z - CUSTOMER SUPPORT

**FOR ANY REQUEST FOR
TECHNICAL SERVICE FOR THE PRODUCT,
WHETHER IN OR OUT OF WARRANTY,
DIRECTLY CONTACT THE

TECHNICAL SUPPORT DEPARTMENT

OF THE DEALER OR RESELLER
THAT SUPPLIED THE PRODUCT.**

ENGLISH

M.O.COM. Srl is completely available to customers to provide any technical information about the product as well as to offer suggestions and advice on steam sterilization procedures.

In this regard, please refer to the following address:

M.O.COM. Srl
Customer Support
Via delle Azalee, 1
20090 Buccinasco (MI)
ITALY
Tel. (+39) 02-45701505
Fax (+39) 02-45701258
e-mail at@mocom.it
website www.mocom.it

To help us in the indispensable work of improving the quality of our products and service, please send your comments and/or suggestions by **letter** or **fax** or to the following **e-mail** address:

uc@mocom.it (Marketing/sale department)

Thank you in advance for your valuable assistance.

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